



ophthalmoscopes & otoscopes



ophthalmoscopes & otoscopes	297 - 324
ophthalmoscopes	298 - 300, 309, 316, 322
otoscopes	300 - 301, 310, 316, 322
ophthalmoscope & otoscope sets	302 - 306, 311 - 314, 317 - 321, 323
wall units	305 - 306, 314, 320 - 321, 323
ophthalmoscope & otoscope accessories	306 - 307, 315, 321, 323

ophthalmoscopes & otoscopes

PanOptic Ophthalmoscope

10
YEARS
GUARANTEE



LED
UPGRADE
AVAILABLE
SEE PAGE 307



WelchAllyn®



- PanOptic Ophthalmoscope gives 5 x larger view of the fundus vs standard ophthalmoscopes in an undilated eye
- 25 degree field of view vs standard 5 degree field of view
- Dioptre range: +20 to - 20
- 6 diagnostic beams: small, medium & large spot, slit, cobalt blue & red-free
- Corneal magnifying lens to help detect corneal abrasions
- Patient eyecup maintains a comfortable viewing distance between practitioner and patient, and occludes ambient light
- Compatible with all existing Welch Allyn 3.5v power sources
- OPH274 comes with PanOptic head, lithium ion handle, charger and hard case

Panoptic iExaminer adapter


The new Welch Allyn iExaminer turns the PanOptic Ophthalmoscope into a mobile digital imaging device allowing you to view and take pictures of the eye like you've never done before.

The adapter aligns the optical access of the PanOptic Ophthalmoscope to the visual axis of the iPhone 4 or 4S camera to capture high resolution (5/ 8 megapixel*) pictures of the fundus and retinal nerve.

The iExaminer App (standard) is free of charge to download from the App store and allows you to acquire and store the pictures to a patient file. The Pro version can also be downloaded for £20.99 from the app store and gives you the additional functionality allowing you to email, print (Requires AirPrint enabled printer) and store up to 1,000 patient files.

And, because the PanOptic Ophthalmoscope provides a 5x larger view of the fundus than a standard ophthalmoscope and 25° field-of-view without having to dilate the pupil, the fundus imaging can happen anytime, anywhere.

* Per iPhone Specs

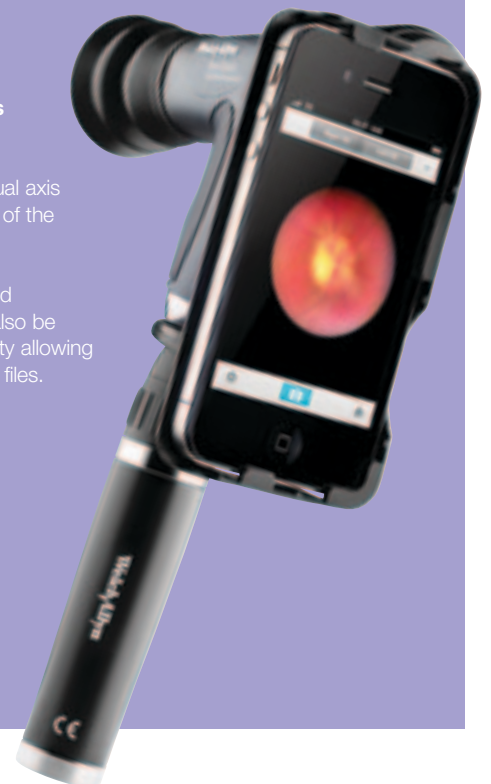

OPH274	Complete Set	£748.00
OPH273	Panoptic iExaminer adapter	£53.00



3.5v Professional Ophthalmoscope

10
YEARS
GUARANTEE



WelchAllyn[®]

- Ophthalmoscope has 6 diagnostic beams: small, medium and large spot, target, slit and red-free
- Dioptre range: +40 to -25
- Comes with c-cell battery handle in a hard case



OPH782		£185.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

3.5v Elite Ophthalmoscope

10
YEARS
GUARANTEE

LED
UPGRADE
AVAILABLE
SEE PAGE 301



WelchAllyn[®]

- Unique coaxial ophthalmoscope allows easy entry to the eye and shadow-free viewing of the retina
- 6 diagnostic beams: small, medium & large spot, target, slit and cobalt blue
- Red-free, and polarising filters
- Dioptre range: +40 to -25
- Double sealed, dust-free optics
- Comes with a c-cell battery handle in a hard case



OPH783		£200.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

ophthalmoscopes & otoscopes

ophthalmoscopes & otoscopes

2.5v Pocketscope Ophthalmoscope



WelchAllyn[®]



- Ophthalmoscope has 6 diagnostic beams: small, medium and large spot, target, slit and cobalt blue, plus red-free filter
- Dioptre range: +40 to -25
- Double sealed, dust-free optics
- Comes with 'AA' metal battery handle and pocket clip in a hard case (batteries not supplied)



OPH474		£167.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04

2.5v Pocketscope Otoscope



- Fibre optic otoscope
- 2.5 x magnification
- Insufflation port for pneumatic otoscopy
- Comes with 'AA' metal handle, pocket clip, throat illuminator, 4 reusable specula and 20 disposable specula in a hard case (batteries not supplied)

WelchAllyn[®]

OPH473		£125.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04



3.5v Professional Otoscope



LED
UPGRADE
AVAILABLE
SEE PAGE 307

10
YEARS
GUARANTEE



- Fibre optic otoscope with HPX bulb for truer tissue colour
- 2.5 x magnification
- Insufflation port for pneumatic otoscopy
- Comes with c-cell battery handle and 4 reusable specula in a hard case

WelchAllyn®

OPH784		£131.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

3.5v Macroview Otoscope



LED
UPGRADE
AVAILABLE
SEE PAGE 307

10
YEARS
GUARANTEE



Better full view of the tympanic membrane than traditional otoscopes - guaranteed



- 4.2 x magnification – 30% greater than a traditional otoscope
- **Twice the field of view** see almost all the tympanic membrane without panning
- **Adjustable focusing** suits different ear canal lengths and compensates for far-sighted eyes
- Insufflation port for pneumatic otoscopy
- Tip grip holds the specula firmly in place and gives ejection mechanism for hands-free disposal
- Comes with c-cell battery handle and 4 reusable specula in a hard case

WelchAllyn®

OPH177		£194.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

ophthalmoscopes & otoscopes

2.5v Pocketscope Set

ophthalmoscopes & otoscopes



WelchAllyn®



10 YEARS GUARANTEE

Ophthalmoscope

- Ophthalmoscope has 6 diagnostic beams: small, medium and large spot, target, slit and cobalt blue, plus red-free filter
- Dioptre range: +40 to -25
- Double sealed, dust-free optics

Otoscope

- Fibre optic otoscope with 2.5 x magnification
- Insufflation port for pneumatic otoscopy
- Comes with 1 x 'AA' metal handle, pocket clip, 4 reusable specula and 20 disposable specula in a hard case (batteries not supplied)

OPH470		£218.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04

3.5v Prestige Diagnostic Set

LED UPGRADE AVAILABLE
SEE PAGE 301

10 YEARS GUARANTEE

WelchAllyn®




Ophthalmoscope

- Unique coaxial ophthalmoscope allows easy entry to the eye and shadow-free viewing of the retina
- Dioptre range: +40 to -25
 - Double sealed, dust-free optics

Otoscope

- Tip grip holds the specula firmly in place and gives ejection mechanism for hands-free disposal
- MacroView otoscope has 4.2 x magnification – 30% greater than a traditional otoscope
- Twice the field of view – see almost all the tympanic membrane without panning
- Adjustable focussing – suits different ear canal lengths and compensates for far-sighted eyes
- Insufflation port for pneumatic otoscopy
- Comes with c-cell battery handle and 4 reusable specula in a hard case

Better full view of the tympanic membrane than traditional otoscopes - guaranteed

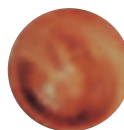

OPH776		£278.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04



3.5v Professional Diagnostic Set



10
YEARS
GUARANTEE

LED
UPGRADE
AVAILABLE
FOR OTOSCOPE

Ophthalmoscope

- 6 diagnostic beams: small, medium and large spot, target, slit and red-free
- Dioptre range: +40 to -25

Otoscope

- Fibre optic otoscope with HPX bulb for truer tissue colour
- 2.5 x magnification
- Insufflation port for pneumatic otoscopy
- Comes with c-cell battery handle and 4 reusable speculum in a hard case



OPH781		£191.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

3.5v Elite Diagnostic Set

LED
UPGRADE
AVAILABLE



Ophthalmoscope

- Unique coaxial ophthalmoscope allows easy entry to the eye and shadow-free viewing of the retina
- 6 diagnostic beams: small, medium & large spot, target, slit and cobalt blue
- Red-free, and polarising filters
- Dioptre range: +40 to -25
- Double sealed, dust-free optics

Otoscope

- Fibre optic otoscope with HPX bulb for truer tissue colour
- 2.5 x magnification and integral throat illuminator
- Insufflation port for pneumatic otoscopy
- Comes with c-cell battery handle and 4 reusable specula in a hard case



10
YEARS
GUARANTEE

OPH675		£230.00
DIS354	Duracell C Size Batteries, pack of 2	£3.04

ophthalmoscopes & otoscopes

3.5v Elite Desk Set

WelchAllyn[®]



**LED
UPGRADE
AVAILABLE**

**10
YEARS
GUARANTEE**

Ophthalmoscope

- Unique coaxial ophthalmoscope allows easy entry to the eye and shadow-free viewing of the retina
- 6 diagnostic beams: small, medium & large spot, target, slit and cobalt blue
- Red-free, and polarising filters
- Dioptre range: +40 to -25
- Double sealed, dust-free optics

Otoscope

- Fibre optic otoscope with HPX bulb for truer tissue colour
- 2.5 x magnification and integral throat illuminator
- Comes complete with either 2 x NICAD or lithium Ion rechargeable handles in a charger



OPH580	3.5v Elite Desk Set	£697.00
OPH581	3.5v Elite Desk Set with Lithium Handle	£944.00

3.5v Prestige Desk Set

WelchAllyn[®]



**10
YEARS
GUARANTEE**

**LED
UPGRADE
AVAILABLE**

Ophthalmoscope

- Unique coaxial ophthalmoscope allows easy entry to the eye and shadow-free viewing of the retina
- 6 diagnostic beams: small, medium & large spot, target, slit and cobalt blue
- Red-free, and polarising filters
- Dioptre range: +40 to -25
- Double sealed, dust-free optics

Otoscope

- 4.2 x magnification – 30% greater than a traditional otoscope
- Twice the field of view – see almost all the tympanic membrane without panning
- Adjustable focusing – suits different ear canal lengths and compensates for far-sighted eyes
- Insufflation port for pneumatic otoscopy
- Comes complete with either 2 x NICAD or lithium Ion rechargeable handles in a charger



OPH778	3.5v Prestige Desk Set	£765.00
OPH789	3.5v Prestige desk set with Lithium Handle	£1033.00

Welch Allyn Green Series™ 777 Wall Unit

A reliable diagnostic solution, configured to your needs, that delivers accurate results whilst providing energy efficiency and cost savings.

- Help reduce costs with a more energy-efficient product than the 767 Wall Transformer
- Help improve workflow and save energy with handles that automatically turn on/off
- Quickly and easily upgrade from the 767 wall transformer to the Green Series 777 wall transformer
- No on/off switch reduces contact, to help reduce the risk of cross-contamination
- 12 foot cables
- Rheostats in handles allow adjustment of light intensity

Standardise or choose the technology you need

Welch Allyn MacroView™ otoscope

- Almost twice the field of view and 30% greater magnification than traditional otoscopes
- Allows for simplified diagnosis of Otitis Media through superior visualisation of the tympanic membrane

Welch Allyn 3.5 V Diagnostic otoscope

- Halogen HPX lamp provides light output for true tissue colour and long-lasting performance
- Fibre optics project cool light to the distal tip with no reflections or obstructions for superior views of the tympanic membrane

Welch Allyn PanOptic™ ophthalmoscope

- Fast, easy entry into undilated pupils with 5X larger view of the fundus for a streamlined exam
- 25° field of view vs. the standard 5° field of view provides optimal visualisation

Welch Allyn Coaxial ophthalmoscope

- 18 unique aperture/filter combinations for greater versatility
- 28 focusing lenses with a range of -25 to +40 diopters



WelchAllyn®

OPH928	GS 777 Wall Transformer	£430.00
OPH930	GS 777 Wall Transformer With Coaxial Ophthalmoscope & Diagnostic Otoscope	£625.00
OPH932	GS 777 Wall Transformer With Coaxial Ophthalmoscope & MacroView Otoloscope	£682.00
OPH933	GS 777 With Panoptic Ophthalmoscope & Macroview Otoloscope	£980.00

Welch Allyn Green Series™ 777 Integrated Wall System

Welch Allyn industry-leading blood pressure and thermometry devices provide you with fast, accurate and reliable patient data you can trust – to help improve patient care and clinical decisions.

All configurations are supplied with the following as standard:

- 777 Wall Transformer
- 767 Aneroid Sphygmomanometer
- Speculum Dispenser
- Easy to install pre-drilled mounting board

Welch Allyn wall aneroid

- Recessed, laser-engraved dial face with high-contrast pointer provides certified accuracy to reduce errors and improve patient care
- Lifetime calibration warranty reduces cost of ownership

Welch Allyn Connex® ProBP™ digital blood pressure device

- Improve patient comfort and increase throughput with fast, 15-second NIBP readings
- Portable, so you can take BP readings just about anywhere
- Connectivity to most leading EMR systems, to improve patient documentation accuracy and efficiency

Welch Allyn SureTemp® Plus electronic thermometer

- Provides accurate oral temperatures in 4 to 6 seconds, axillary and rectal temperatures in 10 to 13 seconds and adult axillary temperatures in 12 to 15 seconds
- Reduce risk of cross-contamination with single-use probe covers
- Improve efficiency with convenient storage for 75 probe covers
- Takes approximately 6,000 readings on 3 AA batteries and provides last temperature recall



WelchAllyn®

OPH937	GS 777 Elite Integrated Diagnostic System Coaxial Ophthalmoscope/Diagnostic Otoloscope	£972.00
OPH935	GS 777 Panoptic Prestige Integrated Diagnostic System Panoptic Ophthalmoscope/ Diagnostic Otoloscope	£1308.00
OPH938	GS 777 Integrated Diagnostic System With Connex PRO BP/Coaxial/Macroview & SureTemp Plus Thermometer	£1776.49
OPH939	GS 777 Integrated Diagnostic System With Connex PRO BP/Panoptic/Macroview & SureTemp Plus Thermometer	£2013.57

ophthalmoscopes & otoscopes

Green Series Elite Diagnostic System

Providing you with accurate, reliable, energy-efficient products you trust in a workflow that will help improve your practice's throughput.

The Green series elite diagnostic system provides the following diagnostic tools:

Rechargeable desk set

- Coaxial Ophthalmoscope
 - 18 unique aperture/filter combinations for greater versatility
 - 28 focusing lenses with a range of -25 to +40 diopters
- Diagnostic Otoscope
 - Halogen HPX lamp provides light output for true tissue colour and long-lasting performance.
 - Fibre optics project cool light to the distal tip with no reflections or obstructions for superior views of the tympanic membrane
- Nicad Rechargeable Handles

767 Aneroid

- Recessed, laser engraved dial face with high-contrast pointer provides certified accuracy to reduce errors and improve patient care
- Lifetime calibration warranty reduces cost of ownership

Pro4000 Thermometer

- ExacTemp™ technology provides an active user feedback system to ensure proper positioning and improved accuracy
- Automatic probe cover eject button for quick and easy disposal of used probe covers



WelchAllyn®

Available
March 2013



OPH676

£826.00

Welch Allyn Spares

WelchAllyn®



Replacement Bulbs

OPH893	Lamp 06500-U for MacroView Otoscope	1	£16.70
OPH871	Lamp 03900 for Pocket Ophthalmoscope		£15.20
OPH872	Lamp 03400 for Clinic Plus & Pocket Otoscope		£15.20
OPH873	Lamp 04900 for 3.5v Elite Ophthalmoscope	2	£18.80
OPH874	Lamp 03100 for 3.5v Elite & Pro Otoscope	3	£16.20
OPH875	Lamp 03000 for 3.5v Halogen Standard Ophthalmoscope	4	£15.20
OPH892	Lamp 03800 for PanOptic	5	£22.00

Otoscope Speculum

OPH878	Set of 4 reusable Specula to fit all Welch Allyn Otoposcopes Sizes 2, 3, 4 & 5mm	6	£12.50
OPH579	Disposable Otoscope Specula. Size 2.75mm or 4.25mm - Specify Size, pack of 850	7	£53.00

Batteries

OPH692	3.5v Nicad Rechargeable Battery, 72300	8	£65.00
OPH877	3.5v Nicad Rechargeable Battery, 72200	9	£61.00
OPH894	3.5v Lithium Ion Rechargeable Battery, 71960	10	£146.00





PanOptic Accessories

OPH278	PanOptic Eyecups, pack of 5		£38.00
--------	-----------------------------	--	--------



Welch Allyn SureColor™ LED Replacement Lamps

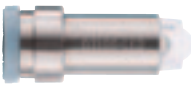7
YEARS
GUARANTEE

LED Lamp Part No.	Device Part No.	Description
03100-LED 	Traditional Diagnostic Oscopes	
	20000	Otoscope/Throat Illuminator
	20200	Pneumatic Oscope
	20200	Pneumatic Consulting Oscope
	20200	Operating Oscope
	20200	Diagnostic Oscope (w/ Specula)
06500-LED 	MacroView™ Diagnostic Oscopes	
	23810	MacroView™ Oscope
	23820	MacroView™ Oscope/Throat Illuminator
	23820	Digital MacroView™ Oscope
04900-LED 	Traditional Ophthalmoscopes	
	11720	3.5V Coaxial Ophthalmoscope
	11730	AutoStep® Coaxial Ophthalmoscope
	11735	Prestige™ Coaxial-Plus Ophthalmoscope
03800-LED 	PanOptic™ Ophthalmoscopes	
	11810	PanOptic™ Ophthalmoscope
	11820	PanOptic™ Ophthalmoscope/Blue Filter/Lens





30X Longer Life

- Tested to 30X halogen lamps
- Guaranteed 7 years



Energy Efficient

- 3X the on-time vs. halogen with Welch Allyn 3.5 V rechargeable handles
- Lower operating temperature intended for greater patient/clinician comfort

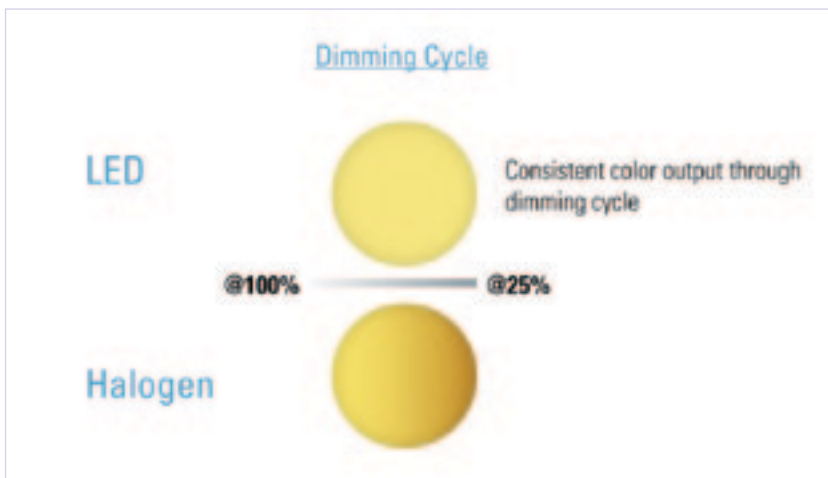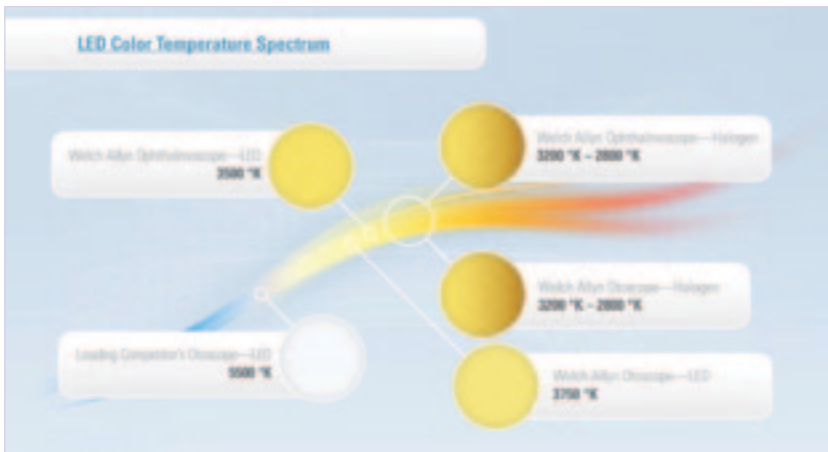

Welch Allyn LED lamps use less energy, which results in lower contact temperature and three times the on-time of your traditional rechargeable handle.

- Patented circuitry that preserves the feel of halogen while offering the longevity and efficient performance of LED technology
- Long life backed by a 7-year warranty ensuring peace of mind and virtually eliminating inconvenient interruptions for lamp replacement
- SureColor™ LED technology enabling reduction of intensity to alleviate glare or enhance patient comfort while preserving the colour and clarity of the view
- High energy efficiency resulting in 3X the examination time with each charge of a Welch Allyn rechargeable handle
- A simple solution presented in four lamp upgrade kits compatible with our most popular physical assessment instruments

OPH881	LED Bulb for Traditional Diagnostic Oscopes	£62.00
OPH882	LED Bulb for MacroView™ Diagnostic Oscopes	£64.00
OPH883	LED Bulb for Traditional Ophthalmoscopes	£72.00
OPH884	LED Bulb for PanOptic™ Ophthalmoscopes	£80.00

Pure brilliance










New Keeler Jazz Diagnostic Set

- Slimline size fits neatly in your pocket
- Trusted Keeler optics
- Brilliant LED Fibre Optic Otoscope
- Competitively priced



Keeler

Apertures, Graticules & Filter Symbols

	<p>Wide Angle Illuminates the largest area of fundus for the best possible general diagnosis through a dilated pupil.</p>		<p>Slit Used primarily to determine retinal elevation and depression, but may also be used to access anterior chamber depth.</p>		<p>Semi Circle Provides a combination of depth perception and field of view.</p>
	<p>Intermediate Permits easier access through an undilated pupil and in peripheral examination. Particularly useful in pediatric examination.</p>		<p>Glaucoma Projects a graticule onto the retina to assess the optic cup/disc ratio as an aid to glaucoma diagnosis.</p>		<p>Red Free The red free filter is used to examine blood vessels in fine detail. By filtering out the red light, blood vessels are silhouetted black against a dark green background.</p>
	<p>Macular Designed specifically for examination of the macular region of the fundus where a large beam would create excessive pupillary reaction or patient discomfort.</p>		<p>Fixation Cross Projects a graticule onto the retina for assessment of the degree and direction of eccentric fixation.</p>		<p>Cobalt Blue Used in conjunction with fluorescein dye for the detection and examination of corneal scars and abrasions.</p>



Practitioner Ophthalmoscope

- The positive maywheel lens control allows you to select from a range of lens powers from +40D to -25D. This increased lens range allows detailed examination of the ocular structures and guarantees ease of use for spectacle wearers.
- The Practitioner includes 6 beams for a full diagnostic capability. A blue filter, particularly useful for assessing corneal scratches in children, is included.
- Red free filter for enhanced vessel examination can be used with the diagnostic beams. This feature is particularly useful for diabetic assessment and monitoring.
- Supplied in hard case with spare bulb
- Requires 2 x AA batteries (not supplied)

OPH267	3.6v Practitioner Ophthalmoscope	£305.00
OPH311	2.8v Practitioner Ophthalmoscope	£215.00
OPH327	Lithium Mini Charger	£115.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH853	2.8v Xenon Ophthalmoscope Bulb, pack of 2	£55.00
OPH854	3.6v Ophthalmoscope Bulb, pack of 2	£55.00



Keeler



2.8v Standard Ophthalmoscope

- 2.8v illumination
- The positive maywheel lens control allows you to click through lens powers from +40D to -25D. This increased lens range allows detailed examination of the ocular structures and provides ease of use for spectacle wearers.
- Two diagnostic beams - a wide angle beam for general examination and a macular beam for more detailed examination of the central retina
- Red free filter for enhanced vessel examination can be used with both diagnostic beams. This feature is particularly useful for diabetic assessment and monitoring.
- The integral soft brow rest allows for easy positioning and protects your lenses if you are a spectacle wearer
- Supplied in hard case with spare bulb
- Requires 2 x AA batteries (not supplied)

OPH310	Standard Ophthalmoscope	£195.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH853	2.8v Xenon Ophthalmoscope Bulb, pack of 2	£55.00



Keeler



2.8v Pocket Ophthalmoscope

- 2.8v illumination
- Full diagnostic capability in a pocket sized instrument
- Comprehensive lens range from +20D to -20D
- Ideal for domiciliary or ward use or as a 'second' instrument
- Security assured from the strong metal pocket clip
- A total of 6 diagnostic beams
- Supplied in hard case with spare bulb
- Requires 2 x AA batteries (not supplied)

OPH250	Pocket Ophthalmoscope	£215.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH850	Pocket Ophthalmoscope Bulb, pack of 2	£55.00



Keeler



ophthalmoscopes & otoscopes

2.8v Pocket Otoscope



Keeler

- 2.8v illumination
- Ideal for domiciliary or ward use or as a second instrument
- Closed head design particularly suited to pneumatic testing
- 2.5x magnifying lens removable to allow the introduction of small instruments
- 4 x Reusable specula supplied
- Security from the strong metal clip
- Supplied in a hard case with spare bulb and bulb extractor
- Requires 2 x AA batteries (not supplied)
- See page 315 for specula

OPH150	Pocket Otoscope	£130.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH855	Pocket Standard Otoscope Bulb, pack of 2	£55.00

2.8v Standard Otoscope



Keeler

- 2.8v illumination
- Closed head design particularly suited to pneumatic testing
- 2.5x magnifying lens removable to allow the introduction of small instruments
- 4 x Reusable specula supplied
- Supplied in a hard case with spare bulb, pneumatic testing adaptor and bulb extractor
- Requires 2 x AA batteries (not supplied)
- See page 315 for specula

OPH305	Standard Otoscope	£135.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH855	Pocket Standard Otoscope Bulb, pack of 2	£55.00

2.8v Practitioner F.O. Otoscope

Keeler



- Fibre Optic 2.8v illumination
- Genuine Halogen illumination transmitted via optic fibres for homogenous illumination leading to clear and accurate diagnosis for aural images
- Large head design for totally unobstructed view of aural structures.
- Closed head system ideally suited to pneumatic testing for tympanic mobility
- Large slide aside 2.5x magnifying lens for clear view of aural structures and wax and foreign body removal procedures under magnification
- Disposable specula supplied and 4 x reusable specula supplied
- Supplied in a hard case with spare bulb, pneumatic testing adaptor
- Requires 2 x AA batteries (not supplied)
- See page 315 for specula

OPH301	Practitioner F.O. Otoscope	£190.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH857	2.8v F.O. /Practitioner Otoscope Bulb, pack of 2	£55.00

2.8v Practitioner Otoscope



Keeler

- 2.8v illumination
- Unobstructed open style head for easy introduction of small instruments for wax or foreign body removal
- Large slide aside 2.5x lens for clear view of aural structures and procedures under magnification
- Facility to undertake pneumatic testing for tympanic mobility
- Disposable specula and 4 x reusable specula supplied
- Supplied in a hard case with spare bulb and airtight chamber/ pneumatic testing adaptor
- Requires 2 x AA batteries (not supplied)
- See page 315 for specula

OPH302	Practitioner Otoscope	£160.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04
OPH857	2.8v F.O. /Practitioner Otoscope Bulb, pack of 2	£55.00

2.8v Pocket Set

Keeler

Pocket Ophthalmoscope

- Full diagnostic capability in a pocket sized instrument
- Ideal for domiciliary or ward use or as a second instrument
- Security assured from the strong metal pocket clip
- A total of 6 diagnostic beams
- Pocket instruments are powered by standard AA size batteries



Pocket Otoscope

- Ideal for domiciliary or ward use or as a second instrument
- Closed head design particularly suited to pneumatic testing
- 2.5x magnifying lens removable to allow the introduction of small instruments
- Security from the strong metal clip

The set is supplied in a hard case with 1 x slimline handle, (batteries not supplied), 1 x replacement ophthalmoscope bulb, 1 x replacement otoscope bulb, 1 x bulb extractor for screw in bulb, 4 x otoscope reusable specula and 1 x otoscope nasal specula.



Ideal for the Practitioner on the move

OPH453	Pocket Set	£230.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04

Replacement Bulbs - Pack of 2

OPH850	Pocket Ophthalmoscope Bulb	£55.00
OPH855	Pocket Standard Otoscope Bulb	£55.00

Practitioner with Fibre Optic Otoscope Set

Practitioner Ophthalmoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell 2.8v

- The positive maywheel lens control allows you to select from a range of lens powers from +40D to -25D. This increased lens range allows detailed examination of the ocular structures and guarantees ease of use for spectacle wearers
- The Practitioner includes 6 beams for a full diagnostic capability.
- A Blue filter, particularly useful for assessing corneal scratches in children, is included
- A Red free filter for enhanced vessel examination can be used with both diagnostic beams. This feature is particularly useful for diabetic assessment and monitoring
- Soft integral brow rest allows for easy positioning and protects your lenses if you are a spectacle wearer



Fibre Optic Otoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell option 2.8v

- Genuine Halogen illumination transmitted via optic fibres for homogenous illumination leading to clear and accurate diagnosis for aural images
- Large head design for totally unobstructed view of aural structures
- Closed head system ideal for pneumatic testing for tympanic mobility
- Large slide aside 2.5x magnifying lens for clear view of aural structures and wax and foreign body removal procedures under magnification

The set is supplied in a hard case with 1 x slimline handle, (batteries not supplied in the Dry Cell version, Lithium battery supplied with Lithium set), 1 x replacement ophthalmoscope bulb, 1 x replacement otoscope bulb, 1 x Pneumatic testing adaptor, 5 x otoscope reusable specula and 4 x otoscope disposable specula.

The Lithium Set is not supplied with charger.



The Primary Care Specialist's choice

Keeler

OPH325	Practitioner F.O. Set, Dry Cell	£330.00
OPH324	Practitioner F.O. Set, 3.6v Lithium, without charger	£420.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04

Replacement Bulbs - Pack of 2

OPH854	3.6v Xenon Ophthalmoscope Bulb	£55.00
OPH858	3.6v F.O./Practitioner Otoscope F.O. Bulb	£55.00
OPH853	2.8v Xenon Ophthalmoscope Bulb	£55.00
OPH857	2.8v F.O./Practitioner Optic Bulb	£55.00
OPH327	Lithium Mini Charger	£115.00
OPH328	Double Lithium Charger/Slim Line	£240.00

ophthalmoscopes & otoscopes

Standard Set

Standard Ophthalmoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell 2.8v

- The positive maywheel lens control allows you to click through lens powers from +40D to -25D. This increased lens range allows detailed examination of the ocular structures and provides ease of use for spectacle wearers
- Two diagnostic beams - a wide angle beam for general examination and a macular beam for more detailed examination of the central retina
- Red free filter for enhanced vessel examination can be used with both diagnostic beams. This feature is particularly useful for diabetic assessment and monitoring
- The integral soft brow rest allows for easy positioning and protects your lenses if you are a spectacle wearer



The Lithium Set is not supplied with Charger.

Standard Otoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell 2.8v

- Closed head design particularly suited to pneumatic testing
- 2.5x magnifying lens removable to allow the introduction of small instruments
- Disposable specula supplied

The set is supplied in a hard case with 1 x slimline handle, (batteries not supplied in the Dry Cell version, Lithium battery supplied with Lithium Set), 1 x replacement ophthalmoscope bulb, 1 x replacement otoscope bulb, 1 x Pneumatic testing adaptor, 1 x bulb extractor for screw-in bulbs, 5 x otoscope reuseable specula.

OPH320	Standard Set, Dry Cell	£255.00
OPH321	Standard Set, 3.6v Lithium, without charger	£340.00
OPH326	Standard Fibre Optic Otoscope Diagnostic Set 2.8V	£235.00
DIS352	Duracell AA size Battery, pack of 4	£3.04
OPH327	Lithium Mini Charger	£115.00
OPH328	Double Lithium Charger	£240.00

Replacement Bulbs - Pack of 2

OPH854	3.6v Xenon Ophthalmoscope Bulb	£55.00
OPH856	3.6v Standard Otoscope Bulb	£55.00
OPH853	2.8v Xenon Ophthalmoscope Bulb	£55.00
OPH855	2.8v Pocket/Standard Otoscope Bulb	£55.00

Colour Coded Grips for Diagnostic Sets



- Easy slip on coloured grip
- Fits all Keeler handles
- Single grip

OPH33401	Pink	£8.00
OPH33403	Blue	£8.00

Practitioner Set

Practitioner Ophthalmoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell 2.8v

- The positive maywheel lens control allows you to select from a range of lens powers from +40D to -25D. This increased lens range allows detailed examination of the ocular structures and guarantees ease of use for spectacle wearers
- Includes 6 beams for a full diagnostic capability. A blue filter, particularly useful for assessing corneal scratches in children, is included
- Red free filter for enhanced vessel examination can be used with the diagnostic beams. This feature is particularly useful for diabetic assessment and monitoring

OPH322	Practitioner Set, Dry Cell	£300.00
OPH323	Practitioner Set, 3.6v Lithium, without charger	£390.00
DIS352	Duracell AA size Battery, pack of 4	£3.04
OPH327	Lithium Mini Charger	£115.00
OPH328	Double Lithium Charger	£240.00

Replacement Bulbs - Pack of 2

OPH854	3.6v Xenon Ophthalmoscope Bulb	£55.00
OPH858	3.6v F.O. /Practitioner Otoscope Bulb	£55.00
OPH853	2.8v Ophthalmoscope Bulb	£55.00
OPH857	2.8v F.O. /Practitioner Otoscope Bulb	£55.00

Practitioner Otoscope - Slim Line handle provides choice of Lithium 3.6v or Dry Cell 2.8v

- Unobstructed open style head for easy introduction of small instruments for wax or foreign body removal
- Large slide aside 2.5x lens for clear view of aural structures and procedures under magnification
- Facility to undertake pneumatic testing for tympanic mobility
- Disposable specula supplied

The set is supplied in a hard case with 1 x slimline handle, (batteries not supplied in the Dry Cell version, Lithium battery supplied with Lithium Set), 1 x replacement ophthalmoscope bulb, 1 x replacement otoscope bulb, 1 x Airtight chamber/Pneumatic testing adaptor, 4 x otoscope disposable specula and 5 x otoscope reusable specula.

The Lithium Set is not supplied with Charger.

Keeler





Jazz Diagnostic Set

Jazz Otoscope

- Fibre Optic Otoscope
- Brilliant Xenon illumination
- 3 x magnifying lens

Jazz Ophthalmoscope

- Brilliant Xenon illumination
- +/- 20D range of powers
- Renowned high performance Keeler optics with aspherical lenses

Jazz Diagnostic Set

- Jazz instruments in one handy package with two handles
- All Jazz instruments are supplied with 3 interchangeable coloured rings



OPH306	Jazz Otoscope	£88.00
OPH312	Jazz Ophthalmoscope	£125.00
OPH335	Jazz Diagnostic Set	£205.00



3.6v Rechargeable Set

1.5% discount
when payment received with order



Desk Unit Sets - Lithium Charger & Slim Line handles

- Ideal for ward, surgery, consulting room and casualty use
- 3.6v Illumination
- Available with Keeler range of ophthalmoscopes and otoscopes
- Can be desk or wall mounted to suit your particular requirements
- Lithium-ion technology provides **over 2 hour use time and 2 hour fast charge**
- LEDs indicate when battery life is low, charging or fully charged
- Complete with two handles and charger
- Choice of:

- Set A** Standard Ophthalmoscope, Standard Otoscope
- Set B** Practitioner Ophthalmoscope, Practitioner Otoscope
- Set C** Practitioner Ophthalmoscope, Practitioner Fibre Optic Otoscope

OPH350	Rechargeable Set A Standard Set	1	£675.00
OPH351	Rechargeable Set B Practitioner Set		£715.00
OPH352	Rechargeable Set C Practitioner with Fibre Optic Otoscope Set		£805.00
OPH329	Charger Wall Mounting Kit		£27.00




ophthalmoscopes & otoscopes

3.6v Cord Handle Set

- Keeler's anti-theft designed ophthalmoscope and otoscopes have a lockable head and handle fitting
- 3.6v Illumination
- The new, compact modern wall mounted design is ideal for ward, surgery, consulting room and casualty
- Strong coiled tangle free cable permits instrument use up to 3 metres from transformer position
- Automatic illumination when the instruments are removed from their cradle. With LED confirmation
- Optional wall mounted Dispos-A-Spec, specular dispenser. Supplied with 25 x 2.5mm, 3.5mm, 4.5mm and 5.5mm disposable specula sizes
- Mains only

Keeler

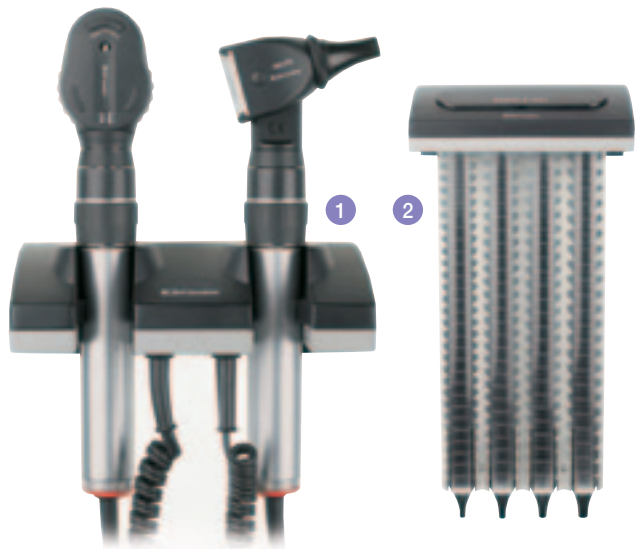

Keep your Instruments where you left them!

OPH353	Cord Handle Set A Standard Set		£695.00
OPH354	Cord Handle Set B Practitioner Set		£720.00
OPH355	Cord Handle Set C Practitioner with F.O. Otoscope set	1	£815.00
OPH340	Dispos-A-Spec	2	£65.00

Windsor BPM Diagnostic Wall Unit



Compact, wall mounted design

- Large Diameter white dial
- Supplied with adult velcro cuff, latex free bulb
- Cuff stored behind dial
- Range 0-300mmHg

OPH356

£950.00




Keeler Bulbs

- All Keeler instruments use genuine Halogen and Xenon technology, yielding a brighter, whiter illumination and a longer life
- Look at the bottom of your bulb to identify if your instrument is 2.8v dry cell battery (identified by a blue circle) or 3.6v rechargeable battery (identified by a red circle)
- Supplied in packs of 2

Ophthalmoscope Bulbs

OPH850	Pocket Ophthalmoscope Bulb 2.8v Halogen	1	£55.00
OPH853	New Standard, New Practitioner Ophthalmoscope 2.8v Bulbs	2	£55.00
OPH854	New Standard, New Practitioner Ophthalmoscope 3.6v Bulbs		£55.00
OPH851	Vista & Vista 20 Ophthalmoscope Bulb 2.8v Halogen	3	£55.00
OPH852	Vista & Vista 20 Ophthalmoscope Bulb 3.6v Halogen		£55.00

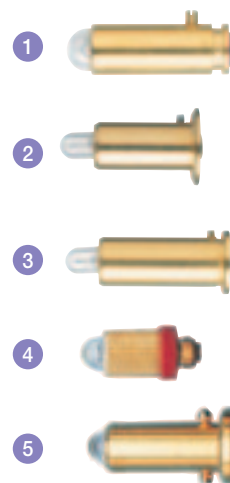

Otoscope Bulbs

OPH855	Standard, Pocket, Deluxe, Vetscope Otoscope & Tongue Depressor Bulb - 2.8v Halogen		£55.00
OPH856	Standard, Deluxe, Vetscope Otoscope and Tongue Depressor Bulb - 3.6v Halogen	4	£55.00

Practitioner, Vista, Fibre Optic etc. Bulbs - push in

OPH857	Practitioner, Vista, Fibre Optic Otoscope & Finhoff Transilluminator Bulb, 2.8v Halogen	5	£55.00
OPH858	Practitioner, Vista, Fibre Optic Otoscope and Finhoff Transilluminator Bulb, 3.6v Halogen		£55.00

Otoscope Specula


Disposable Specula

- Suitable for Practitioner, Standard, Pocket Vista, Fibre Optic. No adaptor required
- Supplied in packs of 500



Reusable Specula


- Suitable for Practitioner, Vista, Fibre Optic, Standard and Pocket. No adaptor required
- Supplied in packs of 5

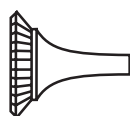

OPH341	2.5mm	£53.00
OPH948	3.5mm	£53.00
OPH342	4.5mm	£53.00
OPH949	5.5mm	£53.00
OPH956	Pack of 12	£9.00

OPH952	Reusable Specula	£27.00
--------	------------------	--------

General Accessories

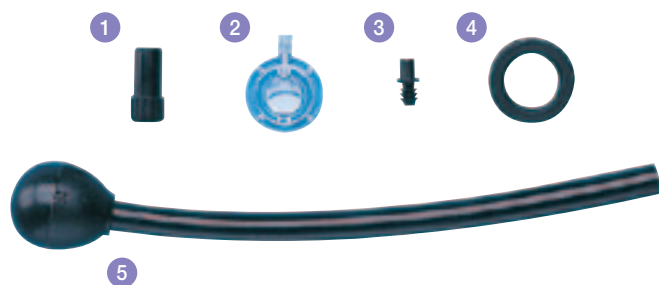

OPH957	Bulb Extractor for Both Old & New Type Standard & Pocket	1	£5.00
OPH958	Insufflator Adaptor for Vista/Practitioner	2	£10.00
OPH960	Insufflator Adaptor for Fibre Optic, New Standard & Pocket	3	£7.00
OPH961	Airtight Chamber/Magnifier for Both Old & New Type Standard & Pocket, round	4	£18.00
OPH165	Otoscope Viewing Lens, Square		£17.00
OPH962	Pneumatic Testing/Insufflator Bulb for all Keeler Oscopes	5	£19.00

Lithium Chargers



- Fast charges in 2 hours
- Provides over 2 hours use time
- LEDs indicate when battery life is low, charging or fully charged

OPH327	Lithium Mini Charger for Slim Line Handle	1	£115.00
--------	---	---	---------



- Recharges two Slim Line handles
- Fast charges in 2 hours
- Provides over 2 hours use time
- LEDs indicate when battery life is low, charging or fully charged
- Handles not included

OPH328	Lithium Charger for Slim Line handle	2	£240.00
OPH329	Charger Wall mounting kit		£27.00

3.75v Lithium Battery



- Lithium battery for use with all new Keeler slim line instrument handles
- When fully charged provides over 2 hours use time

OPH330	Rechargeable Lithium Battery for Slim Line handle	£83.00
OPH331	Red End Cap for Slim Line handle	£26.00

ophthalmoscopes & otoscopes

ophthalmoscopes & otoscopes

2.5v M3000 Ophthalmoscope



- 2.5v illumination
- Lens range +20D to -20D
- Apertures: Large spot, small spot, fixation star, 1/2 circle, red free
- With 18 lenses +/- 1 to -4/+4 in single steps
- High quality handle chrome finished upper section



- On/Off button on clip, auto switch off when returned to pocket
- Supplied in hard case
- Requires 2 x AA batteries (not supplied)

10 YEARS GUARANTEE

OPH204	Ophthalmoscope	£142.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04

2.5v M3000 Otoscopes



- 2.5v illumination
- Choice of fibre optic or direct illumination
- Very bright Xenon illumination
- Swivelling viewing window with 3x magnification
- High quality handle chrome finished upper section
- On/Off button on clip, auto switch off when returned to pocket
- Fibre optic model has insufflation port, (not available on direct illumination model)
- Supplied in hard case with 1 set of re-useable specula, 5 x 2.5mm and 4mm disposable specula
- Requires 2 x AA batteries (not supplied)

10 YEARS GUARANTEE

OPH124	Fibre Optic Otoscope	1	£135.00
OPH125	Direct Illumination Otoscope	2	£87.00
DIS352	Duracell AA Size Batteries, pack of 4		£3.04

2.5v K100/K180 Otoscope



10 YEARS GUARANTEE

- 2.5v illumination, choice of Direct Illumination OPH118 OR Fibre Optic OPH101
- 3x magnification lens movable for insertion of small instruments only on fibre optic model OPH101. Fixed lens on model OPH118
- Closed design
- Suitable for pneumatic testing
- High output Xenon bulb
- Metal handle
- Variable light control rheostat
- Supplied with 3 x reusable specula.
- Complete with soft canvas pouch
- 2 x C size batteries (not supplied)
- Tough ABS construction provides years of reliable use

- Replacement bulbs:
- For K100, select bulb OPH810
 - For K180, select bulb OPH811

OPH101	Fibre Optic Otoscope - K180	£163.00
OPH118	Direct Illumination Otoscope - K100	£128.00
DIS354	Duracell C Size Battery, pack of 2	£3.04

Heine Replacement Bulbs for items on this page

OPH800	2.5v No. 042 - Mini Set & Alpha Ophthalmoscope Bulb	£16.00
OPH813	2.5v No. 056 - Mini Set D.I. Otoscope Bulb	£16.00

OPH810	K100 037 Otoscope Bulb	£16.00
OPH811	K180 077 Otoscope Bulb	£16.00

2.5v M3000 Set



- 2.5v illumination
- Choice of 1 or 2 handle sets
- Choice of fibre optic, OPH124, or direct illumination, OPH125, otoscopes
- M3000 ophthalmoscope, OPH204
- Supplied in hard case with 1 set of
- Reusable specula, 5 x 2.5mm and 4mm disposable otoscope tips
- Each handle requires 2 x AA batteries (not supplied)



OPH442	Direct Illumination 2 Handle Set	£217.00
OPH443	Fibre Optic Illumination 2 Handle Set	£267.00
OPH445	Fibre Optic Illumination 1 Handle Set	£247.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04



2.5v M3000 ENT Multi Set



- ENT diagnostic set
- Compact and durable
- M3000 chrome and plastic combination handle
- M3000 ophthalmoscope head, OPH204
- M3000 Direct Illumination otoscope head, OPH125
- M3000 straight stem laryngeal mirror
- M3000 throat illuminator and tongue depressor
- 4 reusable tips, 2.4/3.0/4.0/5.0mm diameters
- Supplied in hard case
- Requires 2 x AA batteries (not supplied)



OPH446	Multi Set	£249.00
DIS352	Duracell AA Size Batteries, pack of 4	£3.04



OPH805	2.5v No. 105 – M3000 Fibre Optic Otoscope	£16.00
OPH806	2.5v No. 110 – M3000 Direct Illumination Otoscope	£16.00
OPH809	2.5v No. 106 – M3000 Ophthalmoscope	£19.00



Heine Replacement Bulbs for items on this page

ophthalmoscopes & otoscopes

Beta 200 Fibre Optic Diagnostic Set

ophthalmoscopes & otoscopes



Product may vary slightly from image shown

Beta 200 Fibre Optic Diagnostic Set

- Beta 200 ophthalmoscope head
- Beta 200 fibre optic otoscope head
- Set of 4 autoclavable specula
- 10 disposable speculae
- Choice of 2.5v standard battery or 3.5v rechargeable versions
- Supplied complete in hard shell case
- Xenon Otoscope Bulb

Beta 200 Deluxe Fibre Optic Diagnostic Set

- As A332 set, with 3.5v NT300 desk charging system and two 3.5v rechargeable handles (one for each instrument- no changing heads)
- Plus, Heine multipurpose illuminator with tongue blade holder



OPH403	Beta 200 Fibre Optic Diagnostic Set with one battery handle	£501.90
--------	---	---------



2.5v K180 Fibre Optic Set

10
YEARS
GUARANTEE



Set comprises:

- One metal handle
- Variable light control rheostat
- 2 x C size batteries (not supplied)

K180 Ophthalmoscope

- 2.5v illumination
- Lens range +20 to -35D
- Apertures: large, small, red
- free fixation, slit

K180 Fibre Optic Otoscope

- 2.5v illumination
- 3x magnification lens movable for insertion of small instruments
- Closed design
- Suitable for pneumatic testing
- Supplied with 4 x reusable specula

Replacement bulbs:

- For Ophthalmoscope select bulb OPH801
- For Otoscope select bulb OPH811



OPH401	Heine K180 Fibre Optic Set	£310.80
DIS354	Duracell C Size Battery, pack of 2	£3.04

Replacement Bulbs

OPH801	2.5v No. 084 - Kappa 180 Ophthalmoscope Bulb	£16.00
OPH802	2.5v No. 069 - Beta 200 Ophthalmoscope Bulb	£16.00



ophthalmoscopes & otoscopes

K Series 3.5v Fibre Optic Modular Wall Set



1
Standard
OPH541

2
Optional
Instrument
Bay OPH834
with Tongue
Depressor
OPH122

3
Optional
Sphyg SPH200

1 Modular Wall Unit

- Supplied with 2 instrument bays
- Heine 3.5v K180 Ophthalmoscopes
- Heine 3.5v K180 Fibre Optic Otoscope
- 3 Metre spiral cord handles with magnetic automatic on/off and ceramic rheostats
- Instrument locks as standard to prevent unauthorized removal
- Disposable and reusable tips
- Hygienic easy to clean wall mounted unit with smooth rounded corners and no hidden recesses
- Wall mounting Oscopes specula dispenser



OPH541	K Series 3.5v Fibre Optic Wall Set	1	£654.16
--------	------------------------------------	----------	---------

Replacement Bulbs

OPH833	Heine 3.5v No. 086-Kappa 180 Ophthalmoscope Bulb	£20.00
OPH832	Heine 3.5v No. 078-Kappa 180 F.O Otoscope Bulb	£20.00



Optional Modular Equipment

- 2** • 3rd instrument bay, with recommended additional instrument 3.5v tongue depressor with blade ejector and 100 disposable blades
- 3** • XXL Large dial aneroid Sphygmomanometer

3rd Instrument Bay and Instruments

OPH834	3rd Instrument bay and cord handle	2	£159.00
OPH122	3.5v Tongue depressor for disposable blades	2	£66.35
OPH123	3.5v Tongue depressor for wooden tongue blades		£66.35
OPH810	2.5v No. 037- Tongue blade holder Bulb		£16.00

Please note: Any 3.5v Heine instrument can be used provided it has a bayonet fitting.

XXL Large Dial Wall Mounted Aneroid and Available Accessories

SPH200	XXL Large dial aneroid sphygmo with adult cuff and holder	3	£101.92
SPH221	Child cuff for XXL		£28.42
SPH222	Adult cuff for XXL		£23.73
SPH223	Large Adult cuff for XXL		£36.42
SPH224	Thigh Adult cuff for XXL		£49.61
SPH201	Wall mounted cuff holder for additional cuffs		£21.29
SPH225	Inflation bulb with valve for XXL		£9.96
SPH226	Male connector for XXL		£0.49
SPH227	Female connector for XXL		£1.46



Beta 3.5v Fibre Optic Modular Wall Set



1
Standard
OPH541

2
Optional
Instrument Bay OPH834 with Tongue Depressor OPH122

3
Optional
Sphyg SPH200

10 YEARS GUARANTEE

1 Modular Wall Unit

- Supplied with 2 instrument bays
- Heine 3.5v Beta Ophthalmoscope
- Heine 3.5v Beta 200 Fibre Optic Otoscope
- 3 Metre spiral cord handles with magnetic automatic on/off and ceramic rheostats
- Instrument locks as standard to prevent unauthorized removal
- Disposable and reusable tips
- Hygienic easy to clean wall mounted unit with smooth rounded corners and no hidden recesses
- Wall mounting Oscopes specula dispenser



Optional Modular Equipment

- 2 • 3rd instrument bay, with recommended additional instrument 3.5v tongue depressor with blade ejector and 100 disposable blades
- 3 • XXL Large dial aneroid Sphygmomanometer with choice of 1, 2 or 4 additional cuffs and holders

See page 320 for ordering information

OPH542	Beta 3.5v F.O Wall set	£754.95
--------	------------------------	---------

Replacement Bulbs

OPH831	3.5v No. 070 - Beta 200 Ophthalmoscope Bulb	£20.00
--------	---	--------



Heine Oscope Specula and Accessories



OPH826	Insufflation Bulb all Models	£7.40
TOR018	Disposable Tongue Blades 100's	£8.36

Heine Specula

OPH820	Disposable 100's 2.5mm for K100 & Beta 100 D.I. Sets	£10.00
OPH821	Disposable 100's 2.5mm for All Other Oscopes	£10.00
OPH822	Disposable 100's 4mm for K100 & Beta 100 D.I. Sets	£10.00
OPH823	Disposable 100's 4mm for All Other Oscopes	£10.00
OPH824	Autoclavable Set of 6 for Universal D.I. Sets	£47.00
OPH825	Autoclavable Set of 4 for All Other Models	£14.00

Heine Replacement Bulbs for items on this page



OPH814	2.5v No. 041 - Mini Fibrilux F.O. Oscope Bulb	£16.00
OPH810	2.5v No. 037 - Alpha, Mini 2000 F.O. Oscope Tongue Blade Holder Bulb	£16.00
OPH800	2.5v No. 042 - Mini, 2000 & Alpha Ophthalmoscope Bulb	£16.00
OPH813	2.5v No. 056 - Mini, 2000 D.I. Oscope Bulb	£16.00

ophthalmoscopes & otoscopes

ophthalmoscopes & otoscopes

Ri-Scope 2.5v, 3.5v LED and 2.5vHL L2 Otoscope

- Fibre optics for optimal beaming and transmission of the light
- Bi-directional swivelling, sealing optical glass with 3-fold magnification
- Available with Halogen, Xenon or brilliant powerful natural LED



Disposable Ear Specula for Ri-scope L2 and Pen-scope Oscopes

- An especially hygienic option for the prevention of cross-contamination
- For single use
- Made of Blue plastic



OPH184	2.5v HL Otoscope	£148.00
OPH186	3.5v LED L2 Otoscope	£205.00
OPH482	Ri-Scope L2 2.5v HL Set	£302.00
OPH990	2.5mm, pack of 100	£28.50
OPH792	HL 2.5v Halogen Bulbs	£22.80

Ri-Scope 2.5v, 3.5v LED and 3.5HL L2 Ophthalmoscope



- High performance ophthalmoscope with separate red free and cobalt polarisation filters for each aperture
- Dioptre range of +40 to -35 for detailed diagnosis
- Easy to operate aperture hand wheel with six diagnostic beams, semi circle, small medium and large circle fixation star and split for detailed diagnosis
- Optimised high performance optics with aspherical condenser lens
- Parallel beam path
- Dust proof
- Premium hand held diagnostic instruments
- more efficient diagnosis by natural LED illumination
- Patented electronic switch (rheotronic) with automatic switch off
- Durable li-ion batteries
- Renowned high performance/reflection reduced optics (optical glass, aspherical lense)
- 3 years warranty including LED's



OPH283	2.5v HL L2 Ophthalmoscope	£194.00
OPH285	3.5v LED L2 Ophthalmoscope	£262.00
OPH486	Ri-Scope L2 3.5v LED Set	£388.00
OPH793	HL 2.5v Halogen Bulb	£22.80

E-scope

Otoscope

- Economically priced FO version with 2.5 V xenon (3200 K) or 3.7 V LED illumination (5500 K)
- Fibre optics for optimum focusing and guiding of the light
- Sealing precision optics pivotable to both sides with triple magnification for optimum view
- Suitable for performance of pneumatic tests (supplied without connector and insufflator)
- Available in a choice of black or white finish

Ophthalmoscope

- With 2.7 V vacuum or halogen illumination and with innovative 3.7 V LED illumination for higher-contrast and more efficient diagnosis
- Dioptre disc with 18 corrective lenses (+/- 1, 2, 3, 4, 6, 8, 10, 15, 20 dioptres)
- Easy-to-use aperture wheel with six different apertures (fixation star, large circle, small circle, red free filter, blue filter, semi circle) for detailed diagnosis
- High performance optics with aspherical condenser lens
- Spectacle protection
- Parallel optical path
- Dust-protected


• World first and only LED pocket set

- 50 hours service life by using simple alkaline batteries
- Patented technology for consistent light output over service life of batteries
- Brilliant white LED illumination
- Renowned high performance/reflection
- Reduced optics with aspherical lens

OPH18701	Otoscope 2.5v white in case	£66.96
OPH18702	Otoscope 2.5v black in case	£66.96
OPH18703	Otoscop led 3.7v white in case	£85.00
OPH18704	Otoscop led 3.7v black in case	£85.00
OPH28601	Ophthalmoscope 2.5v white in case	£116.50
OPH28602	Ophthalmoscope 2.5v black in case	£116.50
OPH28603	Ophthalmoscope led 3.7v white in case	£135.00
OPH28604	Ophthalmoscope led 3.7v black in case	£135.00
OPH47601	oto/ophtha 2.5v white in case	£145.00
OPH47602	Oto/ophtha 2.5v black in case	£145.00
OPH47603	Oto/ophtha led 3.7v white in case	£175.00
OPH47604	Oto/ophtha led 3.7v white in case	£175.00



L2 Ophthalmoscope and Otoscope Diagnostic Station

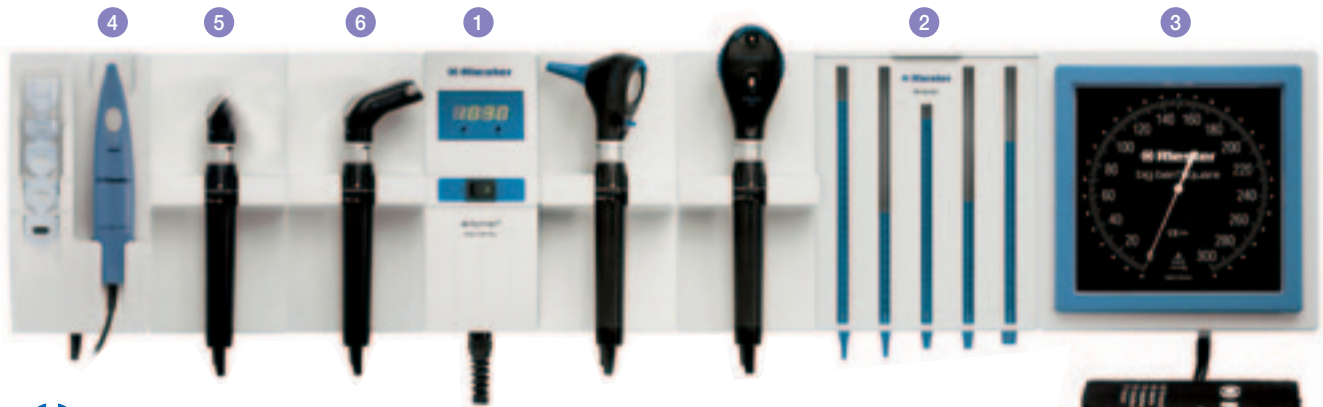

Riester
The familiar way

LED
Lamp

3
YEARS
GUARANTEE
INCLUDING LED

- Anti-theft to ensure your instruments are where you left them
- Can add other module instruments; optional Big Ben sphygmomanometer or Specula dispenser, order separately

- Large scale with 145.5mm long sides. High contrast layout guarantee optimum readability
- Extra sturdy wall mounting made of glass-fibre reinforced plastic and extremely sturdy casing
- Maximum error tolerance: +/- 3mmHg
- Supplied with spacious basket for storage of tubes and cuffs, as well as wall mounting plate, plug and screws

- 3.5v LED lamp
- High quality fibre optic wave conductor for optimum light beam and transmission of light
- Cool examination light similar to daylight
- Slide attachment for fixing and detaching wooden tongue depressor
- Supplied with head only

OPH610	L2 Ophthalmoscope & Otoscope Diagnostic Station	1	£627.00
--------	---	---	---------

SPH170	Big Ben Sphygmomanometer with Adult Cuff	3	£68.50
--------	--	---	--------

OPH612	Tounge Blade	5	£114.00
--------	--------------	---	---------

- Specula dispenser with 5 specula tubes
- Hygienic and practical solution, dispenses a new speculum for each examination
- Simple to clean and disinfect
- Easy to refill with cover, with visible filling levels
- Can also be used independent of the diagnostic station
- Complete with plugs and screws

- Measurement takes approx. 1-2 seconds
- Digital display for optimum readability of the measurement values in °C or °F
- Memory function for 12 measurements
- Supplied with 2 x AA batteries for approx. 1000 measurements
- Supplied with 25 probes

- 3.5v LED lamp
- High quality fibre optic wave conductor for optimum light beam and transmission of light
- Cool examination light similar to daylight
- Swivel and removeable lens with 25x magnification
- Controlled expansion of the speculum by turning on the screw adjuster
- Supplied with head only

OPH611	Ear Specula Dispenser	2	£40.00
--------	-----------------------	---	--------

TUN136	Thermometer	4	£102.60
TUN137	Probe covers, 100pk		£22.80

OPH613	Nasal Speculum	6	£85.50
--------	----------------	---	--------

Ri-charger L



- For your rechargeable 3.5V C-handles
- Table or wall mobile
- Practical storage container for two rechargeable lamps and ear specula
- Charging status displays on the device, Yellow: charging in progress Green: charging completed
- For all rechargeable handles with Li-ion batteries, 120V + 230V

OPH922	ri-charger® L Charging Base	£171.00
--------	-----------------------------	---------

Riester Accessories

For ri-scope F.O. otoscope L2/L3, re-derma, Light conductor, nasal speculum and tongue blade holder



OPH790	2.7v Vacuum Lamp	£22.80
--------	------------------	--------

For Pen-scope Otoscope



OPH791		£22.80
--------	--	--------

LED 3.5V. For ri-scope® F.O. otoscope L2/L3, ri-derma®, Light conductor, nasal speculum and tongue blade holder



OPH993	3.5v Vacuum Lamp	£22.80
--------	------------------	--------

For ri-scope® L1, L2, L3 ophthalmoscope



OPH994	XL 3.5v	£22.80
--------	---------	--------

Disposable Ri-Mini otoscope specula



OPH99201	4mm 500 pack	£55.00
OPH99202	2.5mm 500 pack	£55.00

save time
and money
when you
shop online



 freephone

0800 0304169

 email

info@henryscheinmedical.co.uk

 freefax

0800 0325610

 internet

www.henryscheinmedical.co.uk

 business hours

We are here to take your orders, provide product information and advice. Please ask us about our latest special offers and great deals we would love to offer you.

Monday to Friday
8:30 - 5:30

- The Henry Schein Medical website offers enhanced searching capabilities to help you find what you need faster
- We have centralised all your ordering tools in one location found within 'supplies'
- Pay for orders direct onsite using secure payment technology
- Print copy invoices online and track orders, including proof of delivery
- Request returns and repairs online
- View our current special offers, supplier promotions and clearance lines
- Additionally, all your specific account information will be found within 'my account'
- If you're a new Henry Schein Medical customer or do not already have a web account, it's easy to create a new account simply go to 'create an online account' to set one up